




**A+A Diagnostics**

**Innovative technology for your liquid handling applications**

# BGS - 170

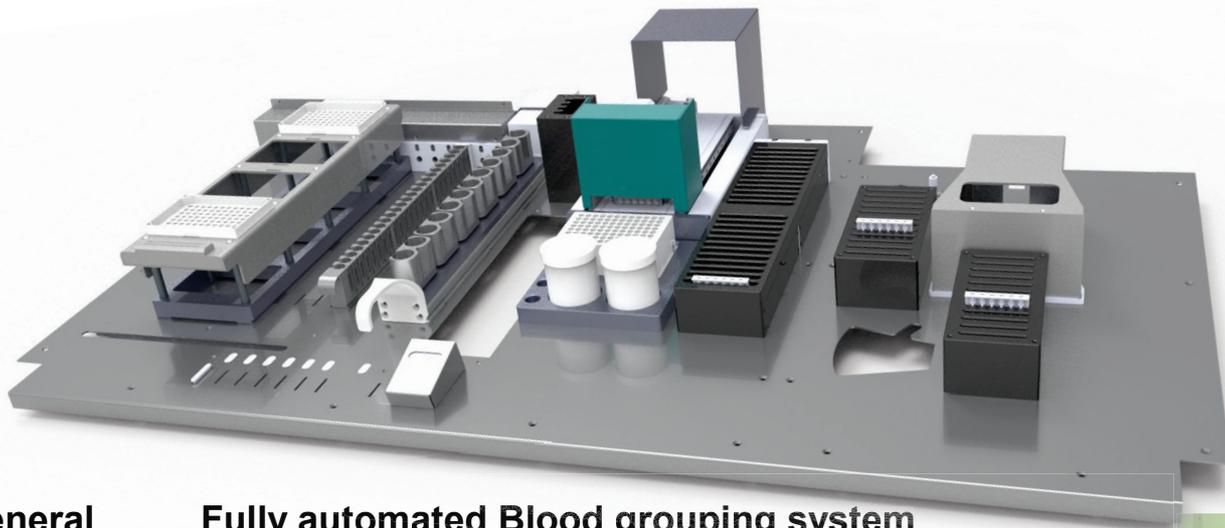

## General Fully automated Blood grouping system

<b>Test Menu</b>	Blood grouping/ reverse grouping/ phenotyping/ crossmatching/ antibody screening/ antibody identification/ direct AHG test(DAT)/ single antigens/Rh(D) Rh(CcEe)/Newborn-ABO-Rh/ Rh/k/ ABO DD/Kell/ AHG Polyspecific/ ADK
<b>Typical</b>	up to 96 cards/ hour
<b>Transfer</b>	1 (transfer card)
<b>Pipetting channels</b>	2 (800ul disposable plastic tips, non conductive, low cost, Liquid level, clot, empty tube detection)
<b>Piercing module</b>	12 cards position
<b>Tips rack capacity</b>	384 tips (4 positions, each can hold 96tips)
<b>Pipetting range</b>	Range: 10 to 800ul/ Increment: 1ul
<b>Sample tube</b>	24 X 6 = 144 samples (standard tube diameter 13mm, other diameters on
<b>Sample types</b>	centrifuged whole blood/ packed red blood cells/ 3-5% and 0.8% red cell suspensions/serum or plasma
<b>Reagent module</b>	2 X6=12 blending diluent position+96 deepwell dilution microplate+2 dilution position+2control
<b>Incubation module</b>	12X2=24 cards position (RT- 50°C adjustable at an increment of 0.1°C )
<b>Centrifuge</b>	12X2=24 cards position (maximum speed 3000r/min, Centrifugal time, speed adjustable programming)
<b>Gel card box</b>	24X4=96 card capacity
<b>Barcode scanner</b>	2 (1 for sample and reagent, 1 for gel card)
<b>Reader</b>	High resolution image acquisition interpretation system, automatic interpretation of acquisition and recognition results.The background of back light system, original image is clear and permanent, easy to query and tracking.
<b>Software</b>	<p>Completely intellectual property rights, humanization design, Wizard-style user interface, easy to use</p> <p>Bi-directional LIS/HIS connection, evaluate and send results</p> <p>Software interface synchronous operation, synchronous refresh, visual display running</p> <p>Operating process real-time monitoring, early warning function with abnormal situation</p> <p>Perfect function of experiment process tracking function of full trace ability tracking analysis and quality control.</p> <p>Powerful memory storage support testing result recheck include original picture, sample number, reagent message,reading result,operator,testing time,test protocol, test method and etc.</p>
<b>Certification</b>	ISO 13485:2003 ISO9001:2008



**A+A Diagnostics GmbH**

Auf der Heide 15 • D- 37351 Dingelstädt

Federal Republic of Germany

P: +49 (0) 36075-439-342 Fx -308

e-mail: info@aa-diagnostics.de

URL: http://www.aa-diagnostics.de